concept

LED maska na obličej

LED maska na tvár

Maska na twarz LED

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LED sejas maska

LED face mask

LED-Gesichtsmaske

Masque facial LED

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Máscara facial LED

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MO8303

CZ

SK

PL

HU

LV

EN

DE

FR

IT

ES

RO

Poděkování

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Acknowledgements

Thank you for purchasing a Concept product. We wish you much joy with your new appliance every day you use it.

Please read the instruction manual carefully before initial use. To refer to this manual any time you need to, we recommend you to keep it in a safe place. And please pass it to any future owner of the appliance.

Technical specifications	
Voltage	5V/1.3A Type C
Battery capacity	2500 mAh 3,7 V
Full charge time	1.5 hours
Max. operating time at lowest intensity	150 minutes
Wavelength – Red light	633 nm + 830 nm + 1072 nm
Wavelength – Blue light	460 nm + 830 nm + 1072 nm
Wavelength – Yellow light	590 nm + 830 nm + 1072 nm
Wavelength – Near infrared light	830 nm
Wavelength – Deep near infrared light	1072 nm
Number of LEDs	100 pcs
Number of chips on each LED	4
Operating time (1 cycle)	10 minutes
Material	liquid medical silicone

IMPORTANT SAFETY PRECAUTIONS

- Do not use the device in any other way than as described in this manual.
- Do not use the LED mask if you are pregnant, suffer from epilepsy, are allergic to light, are taking photosensitizing medications (e.g. some types of antibiotics), or have an active skin infection.
- Before using the device, remove your glasses or contact lenses.
- Remove all the covering and marketing materials from the device before the first use.
- Make sure the mains voltage corresponds to the values on the rating plate of the device.
- Do not leave the device unattended when it is switched on or connected to the mains.
- Before connecting or disconnecting the device from the mains socket, make sure that the device is switched off.


- When disconnecting the device from the wall outlet, never pull the power cord, grasp the plug end and disconnect it by removing it from the receptacle.
- Do not allow children or unauthorized people to handle the device, use it out of their reach.
- This appliance can be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children unless they are aged from 8 years and above and supervised. Children under 8 years of age must be kept away from the appliance and its power supply. Children must not play with the appliance.
- Pay extra attention if the device is used near children.
- Do not use the device as a toy.
- **Before cleaning and after use of the device, switch the device off and disconnect it from the mains.**
- Do not let the power cord touch hot surfaces.
- Do not allow the device to hang on the power cord.
- Do not use the device in an environment where there are explosive gases and flammable substances (solvents, varnishes, adhesives etc.).
- Do not use coarse and chemically aggressive substances to clean the device.
- Do not use the device if it is not working properly, has been dropped, damaged or immersed in liquid. Have it checked and repaired by an authorized service centre.
- Do not use the device with a damaged cord or plug and have it repaired immediately by an authorized service centre.
- Do not wrap the power cord around the device body.
- Do not use the device outdoors.
- The device is only suitable for home use and is not intended for commercial use.
- The device is intended for personal use.
- **Do not immerse the power cord, plug or device in water or any other liquid.**
- Do not use any accessories other than those recommended by the manufacturer.

- Do not repair the device yourself. Contact an authorized service centre.
- The LED mask is not intended for use by children.
- Avoid contact with water. The LED mask is not waterproof.
- Do not store the LED mask in extreme temperatures or humid environments. Store in a cool, dry place out of reach of children.
- Keep your eyes closed during use or wear the provided eye protectors to protect your eyes from the LED light. Direct exposure to the LED light may cause discomfort or potential eye damage.
- If you have any questions, concerns, specific health issues, or experience any side effects, please consult your dermatologist or general practitioner, or contact our customer support.

Damages from failure to follow the manufacturer's instructions are not covered by the warranty.

PRODUCT DESCRIPTION

- 1 LED face mask
- 2 USB cable for mask
- 3 Controller
- 4 On/off button and mode change button
- 5 Brightness level adjustment button
- 6 Cloth bag
- 7 Left and right side straps
- 8 Eye protectors
- 9 USB-C charging cable
- 10 Top strap
- 11 USB-C charging port

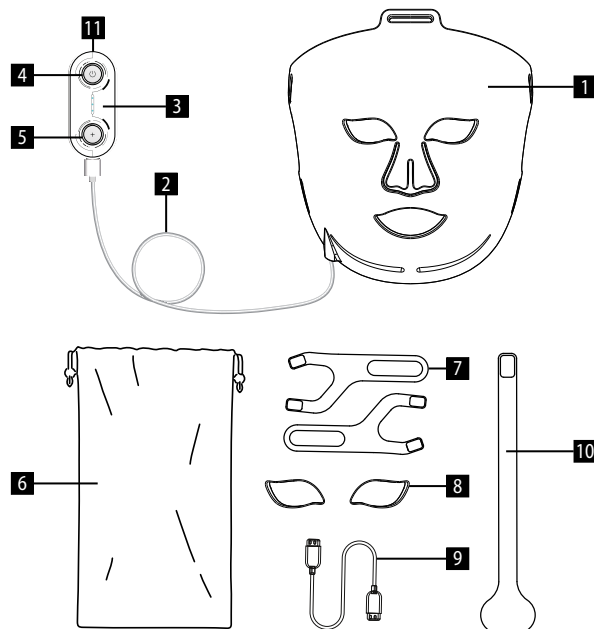INSTRUCTIONS FOR USE

LED light therapy is a biological response of cells to specific wavelengths of light. It is not a laser, but LED radiation, which is safe, does not heat the tissue and is therefore non-invasive.

Some clinical studies published in dermatological literature suggest that light of different wavelengths can promote better skin appearance, such as improving texture or reducing visible imperfections.

Red light (≈620–660 nm) – Penetrates deep into the skin, where it activates fibroblasts, increases collagen and elastin synthesis and can thus reduce fine wrinkles and promote skin elasticity.

Near Infrared Light (≈650–950 nm) – Distinguished by its ability to penetrate deeper into the skin than other visible light wavelengths used in LED therapy. This is a specific range of light wavelengths in the infrared spectrum that help promote deeper tissue regeneration and healing.

Deep near infrared light (≈1000–1042 nm invisible light) – Penetrates deeper layers of tissue (muscles and fascia), improving microcirculation and oxygen uptake by cells. Improves blood circulation, soothes muscle tension and accelerates cellular metabolism for deep skin regeneration and inflammation relief, accelerates healing of acne scars and UV damage.

* Near infrared and deep near infrared light, integrated into LED light therapy can significantly increase the overall effectiveness of therapy.

Blue light (≈405–470 nm) – It acts in the surface layers of the epidermis and, thanks to its antibacterial properties, effectively fights acne-causing bacteria, helps remove imperfections and prevent future acne, reduces inflammation and soothes redness, thereby promoting clearer, smoother skin.

Yellow light (≈570–590 nm) – Penetrates the superficial dermis, the middle layer of the skin, where it helps reduce redness, inflammation, irritation, and leaves the skin refreshed and balanced. Yellow light promotes an even skin tone.

* **Information on the effects of individual wavelengths is based on published scientific studies on LED light in general and may not necessarily correspond to the effects of a specific device.**

Light sensitivity test

- Before first use, perform a light allergy test by placing the LED mask on your forearm for a few minutes. If your skin is red or itchy, do not use the device.
- The LED mask is suitable for all skin types.
- **The skin must be clean, dry and free of makeup or cream before using the LED mask!**

1. Before using the LED mask, thoroughly cleanse your face to remove any dirt, makeup, or cream. This will ensure that the LED light can effectively penetrate your skin. Pat your face dry with a clean towel.
2. Insert the eye shields (8) into the mask and attach the straps to the mask. Thread the ends of the strap (7) through the two holes on the left side and secure the ends in place with the Velcro. Do the same with the second strap on the right side. Continue by threading the strap (10) through the top hole and secure the end in place with the Velcro.
3. Place the LED mask on your face and attach the 3 straps to each other with Velcro.
4. Make sure the controller (3) is charged. Then connect the mask's USB cable (2) to the controller.
5. Using the remote control, turn on the mask by pressing and holding button 4 (approx. 2 seconds).
6. Briefly press button 4 to select the desired mode.
7. Short press button 5 to adjust the brightness level (treatment intensity). The LED mask offers three brightness settings, allowing you to customize the LED light therapy to your skin's sensitivity.
8. Relax during light therapy, keep your eyes closed or wear the included eye shields to protect your eyes from the LED light.
9. If you wish to end the light therapy before the 10-minute timer is up, press and hold button 4 (approx. 2 seconds).
10. There is no need to rinse your face after the LED mask treatment. Apply a moisturizer or serum to your face, which can enhance the effect of the light therapy. LED light therapy can also make the absorption of subsequent skincare products more effective.

Modes


- A Red light 633 nm + 830 nm near infrared light + 1072 nm deep near infrared light** – Skin Enhancement Mode – **Deep near infrared light (1072 nm invisible light)** – Supports the skin's natural processes
- B Blue light 460 nm + 830 nm near infrared light + 1072 nm deep near infrared light** – This mode promotes a clearer appearance of the skin, acting in the surface layers
- C Yellow light 590 nm + 830 nm near infrared light + 1072 nm deep near infrared light** – Mode that helps even out skin tone and promote natural radiance
- * The regimen is suitable for people with sensitive or irritated skin.

Intensity settings

- Low Brightness** – Suitable for those with sensitive skin or those just starting out with LED light therapy. Provides soft and subtle lighting that ensures a pleasant experience without burdening the skin.
- Medium Brightness** – Suitable for people with normal skin sensitivity. Provides an effective dose of LED light that your skin can absorb and benefit from.
- High Brightness** – Suitable for those who prefer stronger stimulation or have thicker skin that may require a higher dose of LED light for effective penetration. High Brightness provides a more intense and brighter LED light output.

Warning: However, it is important to assess your skin's reaction and sensitivity before intensively using high brightness levels, as prolonged exposure to intense light may cause discomfort or irritation.

Warning: It is always a good idea to start with a lower brightness level setting and observe how the skin reacts during and after therapy, and gradually increase the brightness if necessary.

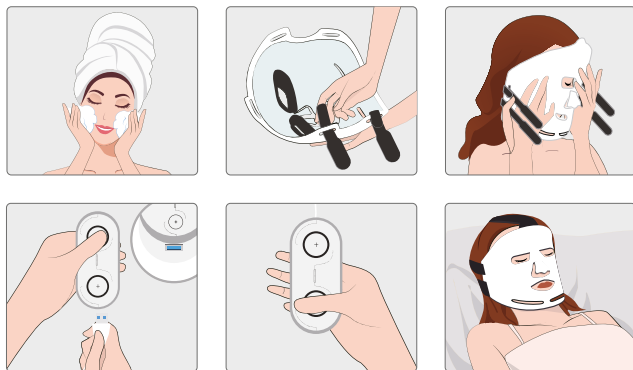Auto-off function - The LED mask has an auto-off function after 10 minutes of operation. This helps prevent overexposure to LED light and ensures that you get the desired benefits from LED therapy and remain safe. For optimal results, we recommend using the LED mask for 10-20 minutes/1 light therapy session. Use regularly 2-4 times a week.

Note: After using LED therapy, your skin may become slightly red, this is normal and should subside within 1-2 hours. If you experience any discomfort, irritation or allergic reactions, please stop using the mask.

Do not exceed the recommended treatment duration and frequency of use to avoid overuse and potential adverse effects on the skin.

Note: Individual results may vary depending on skin type, current condition, and frequency of use.

Charging

- To charge, use the included USB-C charging cable and plug it into the USB-C charging port on the controller. Connect to the adapter and then to a power outlet.
- Recommended adapter parameters: with a maximum output current of 5V/2A.
- Do not wear the LED mask while charging!
- Before charging, we recommend disconnecting the controller from the LED mask and charging it separately.
- Note:** The LED mask cannot be turned on while charging.
- The LEDs on the remote control will light up during charging.
- Full charging time is approximately 1.5 hours, then immediately disconnect from power.

CLEANING AND MAINTENANCE

- After each use, gently clean the inside of the LED mask with a damp soft cloth and then wipe dry. Do not use any cleaning agents or hard objects as they may damage the surface of the appliance!
- Store the LED mask in a dry and dust-free place (ideally in the supplied packaging/fabric bag) away from direct sunlight and extreme temperatures.
- Keep out of reach of children.
- Never clean the appliance under running water, rinse it or immerse it in water or any other liquid!**

TROUBLESHOOTING

Problem	Possible cause	Solution
The controller does not respond while charging.	Damaged charging cable.	Replace with the appropriate charging cable.
		If the problem persists, contact an authorized service center.
The LED on the remote control is flashing.	A flashing LED indicates that the device is discharging.	Charge the device as soon as possible. Let the controller charge until the lights stop flashing and remain solid.
The light in the LED mask does not turn on.	The battery is almost empty.	Fully charge the remote control (approximately 1.5 hours). Note: The LED mask cannot be turned on while charging.
Uneven light output.	Impurities on the surface of LEDs that affect light distribution.	Gently clean the surface with a damp, soft cloth.
	Damaged LEDs.	Contact an authorized service center.
Inconsistent or flickering light.	Check the mask and controller connections.	Charge the controller. If the problem persists, charge the controller in a different outlet or with a different compatible charging cable/adaptor.
Skin irritation, discomfort.	You have exceeded the recommended treatment duration or frequency of use.	Follow the recommended usage instructions – choosing the appropriate light mode, duration, frequency of use, and brightness level.
	Sensitive skin.	Before intensive use of the mask, try it briefly (about 1 minute) with low light intensity.
		If your skin is red or itchy, stop using the LED mask and consult your dermatologist.

SERVICE CENTRE

Extensive maintenance or repair, which requires access to internal parts of the product, has to be performed by a qualified service centre.

ENVIRONMENTAL CONCERNS

- Packaging materials and old appliances should be recycled.
- Packaging materials may be disposed of as sorted waste.
- Dispose the plastic bags made of polyethylene (PE) of as sorted waste.

Recycling of the appliance at the end of its lifespan

This appliance is labelled in accordance with European Directive 2012/19/EU concerning used electrical and electronic appliances (waste electrical and electronic equipment - WEEE). The symbol on the product or on its packaging indicates that this product may not be treated as household waste. Instead it should be taken to the appropriate collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. The disposal of such material has to be done according to the recycling regulations. For more detailed information about recycling of this appliance see the authorized local office or service for household waste disposal or the shop, where you purchased the appliance.

The product meets all the necessary requirements of the EU directives applicable to the specific product.

Changes in text, design, and technical specifications are subject to change without notice and we reserve the right to change them.

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